Digitale Personenwaage Allegra

Allegra digital bathroom scale | Báscula de baño digital Allegra
Pèse-personne numérique Allegra | Bilancia pesapersona digitale
Allegra | Cyfrowa waga łazienkowa Allegra
Allegra digitale personenweegschaal

BE2102

BE2102-1

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BE2102 | BE2102-1

Operating instructions | Instrucciones de uso | Mode d'emploi |
Istruzioni per l'uso | Instrukcja obsługi | Gebruiksaanwijzing

Liebe Kundin, lieber Kunde !

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Dear Customer!

You have chosen to purchase a high-quality ADE brand product, which combines intelligent functions with an exceptional design. This bathroom scale makes it possible for you to determine and monitor your body weight. The wealth of experience of the ADE brand ensures a high technical standard and proven quality.

We hope you enjoy weighing yourself and wish you success in reaching your goal weight!

Your ADE Team

General information

About this operating manual

This operating manual describes how to operate and maintain your product safely. Keep this operating manual in a safe place for future reference. If you pass this product on to someone else, pass on this operating manual as well.

Failure to observe this operating manual may result in injuries or damage to your product.

Explanation of symbols

This symbol combined with the word DANGER warns of serious injuries.

This symbol combined with the word WARNING warns of moderate and mild injuries.

This symbol combined with the word NOTE warns of material damage.

This symbol indicates additional information and general notes.

Table of Contents

General information	10
Safety	11
Intended use	13
Scope of delivery	13
Technical data	13
Overview	14
Start-up	14
Weighing	15
Cleaning	16
Fault/Remedy	16
Declaration of conformity	16
Warranty	17
Disposal	17

Safety

- This product may be used by children aged 8 and above, as well as by persons with reduced physical, sensory, or mental capacities or lacking experience and knowledge if they are supervised or have been instructed in the safe use of the product and understand the risks that may arise.
- Children must not play with the product.
- Cleaning and user maintenance must not be carried out by children, unless they are 8 years of age or older and are being supervised.
- Batteries, when swallowed, can be life-threatening. For this reason, the scale and battery should be kept in a place out of reach of children. If a battery has been swallowed, medical assistance must be sought immediately.

DANGER for children

- Keep packaging material away from children. There is a risk of choking if swallowed. Also pay attention that children do not pull the packaging material over their heads.

RISK of fire/burning and/or explosion

- There is a risk of explosion if the batteries are replaced improperly. Replace the batteries only with the same or equivalent battery type. Pay attention to the "Technical data".
- Batteries must not be recharged, reactivated by other means, dismantled, thrown into fire, or short-circuited. Otherwise, there is an increased risk of explosion and leakage; furthermore, gases may escape!
- Avoid contact of the battery fluid with skin, eyes, and mucosa. In case of contact, immediately rinse the affected areas with plenty of clean water and consult a doctor.

WARNING of risk of injury

- Do not step onto the outer edge of any one side of the scale. Risk of tipping!
- Do not step onto the scale with wet feet or when its surface is damp. Risk of slipping!
- If the scale falls down and shatters, there is a risk of injury from glass splinters.

BEWARE of material damage

- Place the scale on a stable, even surface.
- Do not overload the scale. Pay attention to the "Technical data". Place yourself on the scale exclusively for weighing purposes.
- Protect the scale from strong shocks and vibrations.
- Remove the battery from the scale when it is depleted or when you are no longer using the scale. You will thus avoid damage that may result from leakage of battery acid.
- When inserting and changing the battery, pay attention to the correct polarity (+/-).
- Do not expose the battery to any extreme conditions, for example, by storing it on radiators or under direct sunlight. Increased risk of leakage!
- Clean the battery and device contacts, as needed, before inserting the battery.
- Do not make any modifications to the product. Only allow repairs to be carried out by a specialist workshop. Repairs carried out improperly can considerably endanger the user.
- Do not place the scale in the vicinity of devices that may emit electromagnetic radiation, such as mobile phones, wireless devices, or the like. Doing so may lead to incorrect display indications or malfunctions.

- There are feet on the underside of the scale. Since floor surfaces are made from various materials and are treated with diverse care products, it cannot be completely excluded that some of these substances contain ingredients that may adversely affect and soften the feet of the device. If necessary, place a non-slip pad under the scale.

Intended use

The bathroom scale is designed for the weighing of people in the home. It is not suitable for use in clinics, medical practices, etc.

The maximum carrying capacity of 180 kg must not be exceeded.

Scope of delivery


- Bathroom scale
- 3x AAA / 1.5 V batteries 
- Operating Manual

Technical data

Model:	BE2102, BE2102-1 Allegra digital bathroom scale
Battery:	3x AAA, 1.5 V 
Rated current:	< 40 mA
Measuring range:	max. 180 kg
Variance:	1% +/- 100 g
Display scaling:	100 g
Scale dimensions:	approx. 307 x 307 x 22 mm
Distributor (no service address):	GREND S GmbH Stahltwiete 23 22761 Hamburg, Germany

Our products are constantly being developed and improved.
For this reason, design and technical modifications are possible at any time.

Overview

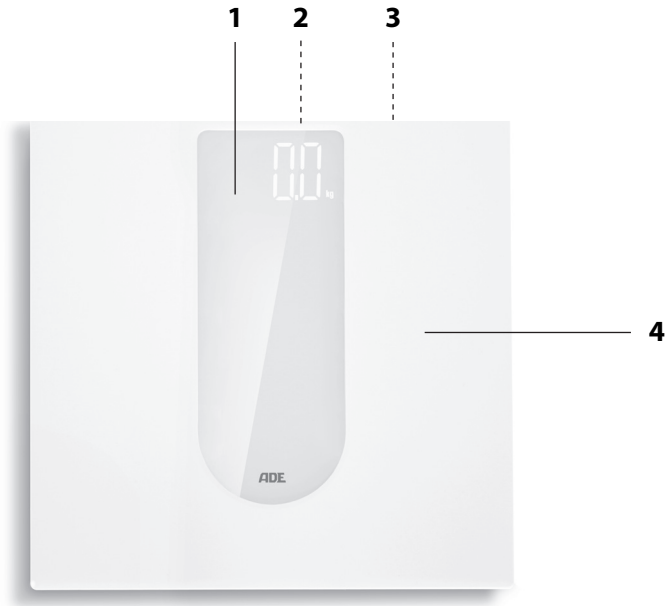- 1 Display
- 2 UNIT button (bottom): Changing the unit of measurement
- 3 Battery compartment (bottom)
- 4 Weighing surface

Start-up

Inserting/changing batteries

For use you need 3 batteries of type AAA, 1.5 V (included). The battery compartment is located on the underside of the device.

1. Push the locking tab of the battery compartment cover towards the battery compartment and remove the cover.
2. Insert the batteries as indicated at the bottom of the battery compartment. Pay attention to the correct polarity (+/-).
3. Replace the battery compartment cover.
The locking tab must snap audibly into place.

Suggestions for changing the battery:

- A change of batteries becomes necessary when “Lo” appears on the display, or when nothing at all is displayed.
- Only use the type of battery that is specified in the “Technical data”.
- Always change all batteries at the same time.
- Dispose of old batteries in an environmentally appropriate way (see “Disposal”).
- Calibrate the scale, see “Setting up the scale”.

Choosing the unit of measurement

The standard unit of measurement upon delivery is “kg” (kilogramme). However, the measurement can also be displayed in the English unit of measurement “lb” (pounds).

1. To turn on the scale, briefly step on the weighing surface or press it.
2. Use the UNIT button on the back of the scale to set the desired unit of measurement. Note the display.

Setting up the scale

1. Place the scale on an even, firm and dry surface. Do not place the scale on a carpet.
2. First, briefly step onto the middle of the weighing surface to calibrate the scale. As soon as the display shows “0.0”, calibration is complete.
This calibration must be performed after each battery change.

Weighing

WARNING against injuries

- Do not step onto the scale with wet feet or when its surface is damp.
Risk of slipping!
- Position yourself on the scale. Stand still and do not hold on to anything.
The measured weight will be displayed shortly thereafter.
A few seconds after you step off the weighing surface, the scale will switch off automatically.

- Weigh yourself with as little clothing on as possible.
- Try to always weigh yourself at the same time and under consistent conditions in order to obtain comparable values.

Cleaning

NOTE on material damage

- Do not immerse the scale in water or other liquids.
- Never use harsh, scouring or abrasive cleaning agents or the like to clean the scale. Otherwise, the surface could be scratched.

1. Wipe the scale as needed with a slightly damp soft cloth.
2. Disinfect the weighing surface from time to time using a suitable disinfecting agent (available at specialist retailers and pharmacies).

Fault/Remedy

Fault	Possible cause
No function	▪ Are the batteries empty or inserted the wrong way around?
The scale appears to be displaying wrong weighing results	▪ Is the scale located on a firm, even surface? ▪ Is a device that can cause electromagnetic disturbances located in the vicinity of the scale? ▪ Remove the batteries from the scale and wait a moment. Reinsert the batteries.
“Err” is shown on the display	▪ The maximum capacity of the scale (180 kg) has been exceeded.
“Lo” is shown on the display	▪ The batteries are empty (see “Inserting/ changing batteries”).

Declaration of conformity

GREND S GmbH hereby declares that the BE BE2102-1 digital bathroom scale complies with the 2014/30/EU and 2011/65/EU directives.

This declaration will lose its validity if any modification is made to the scale without our consent.

Hamburg, june 2025
GREND S GmbH

Warranty

GREND'S GmbH warrants for 2 years from the purchase date the cost-free rectification of defects in material and workmanship through repair or exchange.

In the event of a warranty claim, please return the product with the proof of purchase (specifying the reason for claim) to your dealer.

Disposal

Disposing of the packaging

Dispose of the packaging according to type. Dispose of cardboard with waste paper and of films with recyclable materials.

Disposing of the product

Dispose of the product in accordance with the regulations that apply in your country.

Devices must not be disposed of in normal household waste.

At the end of its life, the product must be disposed of in an appropriate way. In this way, valuable materials contained in the device will be recycled and the burden on the environment avoided. Bring the old device to a collection point for electronic waste or to a recycling depot. For more information, contact your local waste disposal firm or your local administrative body.

Batteries and rechargeable batteries should not be disposed of in the household waste.

All batteries and rechargeable batteries must be brought to a local or district collection point or to the retailer. Such batteries and rechargeable batteries can thus be disposed of in an environmentally friendly way.

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GRENDS GmbH