Personenwaage Olivia

Bathroom Scales Olivia | Báscula personal Olivia
Pèse-personne Olivia | Bilancia pesapersona Olivia |
Waga osobowa Olivia | Personenweegschaal Olivia

Deutsch

English

Español

Français

Italiano

Polski

Nederlands

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BE 1510

Operating Manual | Manual de instrucciones | Mode d'emploi |
Istruzioni per l'uso | Instrukcja obsługi | Gebruiksaanwijzing

Liebe Kundin, lieber Kunde !

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Dear Customer!

You have chosen to purchase a high-quality ADE brand product , which combines intelligent functions with an exceptional design. These bathroom scales make it possible for you to determine and monitor your body weight. The wealth of experience of the ADE brand ensures a high technical standard and proven quality.

We hope you enjoy weighing yourself and wish you success in reaching your goal weight!

Your ADE Team

General information

About this operating manual

This operating manual describes how to operate and maintain your product safely. Keep this operating manual in a safe place for future reference. If you pass this product on to someone else, pass on this operating manual as well.

Failure to observe this operating manual may result in injuries or damage to your product.

Explanation of symbols

This symbol combined with the word DANGER warns of serious injuries.

This symbol combined with the word WARNING warns of moderate and mild injuries.

This symbol combined with the word NOTE warns of material damage.

This symbol indicates additional information and general notes.

Table of Contents

General information	12
Safety	14
Intended use	15
Overview	16
Scope of delivery	16
Starting up	17
Weighing	18
Resetting the scale	18
Cleaning	18
Fault/Remedy	19
Technical data	19
Declaration of Conformity	20
Warranty	20
Disposal	21

Safety

- This product may be used by children aged 8 and above, as well as by persons with reduced physical, sensory, or mental capacities or lacking experience and knowledge if they are supervised or have been instructed in the safe use of the product and understand the risks that may arise.
- Children must not play with the product.
- Cleaning and user maintenance must not be carried out by children, unless they are 8 years of age or older and are being supervised.

Batteries can be life-threatening. Therefore, keep batteries out of reach of small children. If a battery has been swallowed or has entered the body in any other way, seek medical help immediately. Failure to do so within 2 hours may result in severe internal burns that could cause death.

DANGER for children

- Keep packaging material away from children. There is a risk of choking if swallowed. Also pay attention that children do not pull the packaging material over their heads.

RISK of fire/burning and/or explosion

- There is a risk of explosion if the battery is replaced improperly. Replace the battery only with the same or a similar type. Pay attention to the “Technical data”.
- Batteries must not be recharged, reactivated by other means, dismantled, thrown into fire, or short-circuited. Otherwise, there is an increased risk of explosion and leakage; furthermore, gases may escape!
- Avoid contact with skin, eyes, and mucosa. If contact is made with the battery fluid, rinse the affected areas immediately with plenty of clean water and see a doctor straight away.

WARNING of risk of injury

- Do not step onto the outer edge of any one side of the scales. Risk of tipping!
- Do not step onto the scales with wet feet or when its surface is damp. Risk of slipping!

BEWARE of material damage

- Place the scale on a stable, even surface.
- Do not overload the scale. Pay attention to the “Technical data”. Place yourself on the scales exclusively for weighing purposes.
- Protect the scale from strong shocks and vibrations.


- Remove the battery from the scale when it is depleted or when you are no longer using the scale. You will thus avoid damage that may result from leakage of battery acid.
- When inserting and changing batteries, pay attention to the correct polarity (+/-).
- Do not expose the battery to any extreme conditions, for example, by storing it on radiators or under direct sunlight. Increased risk of leakage!
- Clean the battery and device contacts, as needed, before inserting the battery.
- Do not make any modifications to the product. Only allow repairs to be carried out by a specialist workshop. Repairs carried out improperly can considerably endanger the user.
- Do not place the scales in the vicinity of devices that may emit electromagnetic radiation, such as mobile phones, wireless devices, or the like. Doing so may lead to incorrect display indicators or malfunctions.
- There are feet on the underside of the scale. Since surfaces are made from various materials and are treated with diverse care products, the possibility cannot be completely excluded that some of these substances contain constituents that attack and soften the feet of the device. If necessary, place a non-slip pad under the scale.

Intended use

The bathroom scales are designed for the weighing of people in the home. It is not suitable for use in clinics, medical practices, etc.

The maximum carrying capacity of 180 kg must not be exceeded.

Overview

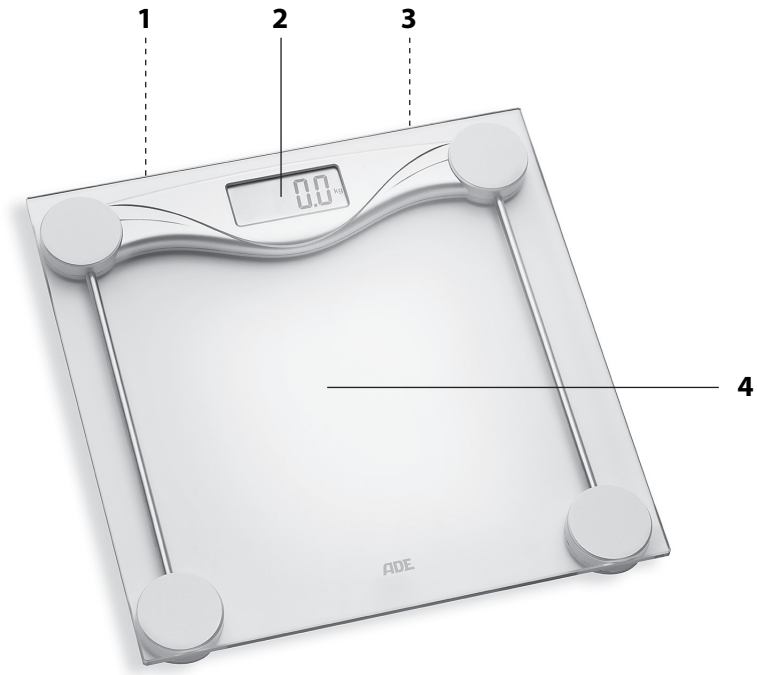- 1 **kg/lb/st** button (underside): Changing the unit of measurement
- 2 Display
- 3 Battery compartment (underside)
- 4 Weighing surface

Scope of delivery

- Bathroom scales
- Battery type CR2032, 3 V 
- Operating Manual

Starting up

Inserting/changing batteries

For use you require one battery of the type CR2032, 3 V (included). The battery compartment is located on the underside of the device.

1. Remove the battery compartment cover.
2. Push the battery towards the contact springs under the big, angular retaining plate and then press it down. The battery will snap into the holding lugs on the opposite side and the battery will be held firmly in its holder.
Pay attention to the correct polarity (+/-). The plus pole (+) must point upwards.
3. Replace the battery compartment cover. The battery compartment cover must snap securely into place.

Suggestions for changing the battery:

- A battery change is necessary when "Lo" or nothing at all is shown on the display.
- Only use the type of battery that is specified in the "Technical data".
- Dispose of the old battery in an environmentally responsible way (see "Disposal").
- Calibrate the scale, see "Setting the scale".

Choosing the unit of measurement

When the scales are delivered, the unit of measurement is set to "kg" as standard. However, the measurement result can also be displayed in the English unit "lb" (pound) or "st" (stone).

- Push the **kg/lb/st** button on the back of the scales.
Note the display.

Positioning the scales

1. Place the scale on an even, firm, and dry surface. Do not place the scales on a carpet.
2. First, briefly step onto the middle of the weighing surface to calibrate the scale. As soon as the display shows "0.0", calibration is complete. This calibration must be performed after each battery change.

Weighing

WARNING against injuries

- Do not step onto the scales with wet feet or when its surface is damp.
Risk of slipping!

- Position yourself on the scales. Remain standing still and do not hold on to anything.

The scale switches on and the display flashes shortly afterwards: The weight measurement is displayed in the selected unit of measure.

A few seconds after you step off the weighing surface, the scales will switch off automatically.

- Weigh yourself with as little clothing on as possible.
- Try to weigh yourself always at the same time and under consistent conditions in order to obtain comparable values.

Resetting the scale

Reset the scale to its default settings if unusual or obviously incorrect values are shown on the display.

1. Remove the batteries from the battery compartment and wait a moment.
2. Re-insert the batteries. Pay attention to the correct polarity of the batteries (+/-). The plus pole (+) must point upwards.

After switching the scales on, you can weigh yourself as usual.

Cleaning

NOTE on material damage

- Do not immerse the scale in water or other liquids.
 - Do not under any circumstances use harsh, scouring or abrasive cleaning agents or the like to clean the scale. Otherwise, the surface could be scratched.
-
- Wipe the scale as needed with a slightly damp soft cloth.
 - Disinfect the weighing surface from time to time using a suitable disinfecting agent (available at specialist retailers and pharmacies).

Fault/Remedy

Fault	Possible cause
No function	Is the battery empty or inserted the wrong way around?
The scale appears to be displaying wrong weighing results	- Is the scale located on a firm, even surface? - Is a device that can cause electromagnetic disturbances located in the vicinity of the scale? - Remove the battery from the scales and wait a moment. Reinsert the battery.
"Err" is shown on the display	The maximum capacity of the scale (180 kg) has been exceeded.
"Lo" is shown on the display	The battery is empty (see "Inserting/ changing the battery").

Technical data

Model:	Bathroom Scales BE 1510 "Olivia"
Battery:	1x CR2032/3 V 
Rated current:	4 mA
Measuring range:	max. 180 kg
Variance:	1 % +/- 100g
Display scaling:	100 g
Weight:	approx. 1.2 kg
Scale dimensions:	approx. 280 x 280 x 24 mm
Distributor: (no service address)	GREND S GmbH Stahlwiete 23 22761 Hamburg, Germany

Our products are constantly being developed and improved.
For this reason, design and technical modifications are possible at any time.

Declaration of Conformity

GREND S GmbH hereby declares that the electronic bathroom scale BE 1510 complies with EU directives 2014/30/EU and 2011/65/EU.

The declaration will lose its validity if any modification is made to the device without our consent.

Hamburg, June 2025
GREND S GmbH

Warranty

GREND S GmbH warrants for 2 years from the purchase date the cost-free rectification of defects in material and workmanship through repair or exchange.

In the event of a warranty claim, please return the product with the proof of purchase (specifying the reason for claim) to your dealer.

Disposal

Disposing of the packaging

Dispose of the packaging with similar materials. Dispose of cardboard with waste paper and of films with recyclable materials.

Disposing of the product

Dispose of the product in accordance with the regulations that apply in your country.

Devices must not be disposed of in normal household waste

At the end of its life, the product must be disposed of in an appropriate way. In this way, valuable materials contained in the device will be recycled and the burden on the environment avoided. Bring the old device to a collection point for electronic waste or to a recycling depot.

For more information, contact your local waste disposal firm or your local administrative body.

Standard and rechargeable batteries should not be disposed of along with household waste

All standard and rechargeable batteries must be brought to a local or district collection point or to the retailer. Such standard and rechargeable batteries can thus be disposed of in an environmentally friendly way.

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GREND'S GmbH