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Dear Customer!

Thank you for purchasing a high-quality ADE brand product that combines intelligent functions with an exceptional design. This bathroom scale makes it possible for you to determine and monitor your body weight. The wealth of experience of the ADE brand ensures a high technical standard and proven quality.

We hope you enjoy weighing yourself and wish you success in reaching your goal weight!

Your ADE Team

General information

About this operating manual

This operating manual describes how to operate and maintain your product safely. Keep this operating manual in a safe place for future reference. If you pass this product on to someone else, pass on this operating manual as well.

Failure to observe this operating manual may result in injuries or damage to your product.

Explanation of symbols

This symbol combined with the word DANGER warns of serious injuries.

This symbol combined with the word WARNING warns of moderate and mild injuries.

This symbol combined with the word NOTE warns of material damage.

This symbol indicates additional information and general notes.

Table of Contents

General information	16
Scope of delivery	17
Safety	18
Intended use	19
Overview	20
Scale.....	20
External Display – Front.....	20
External Display – Back.....	21
Start-up	22
Set up/hang external display.....	22
Using the scale.....	23
Weighing.....	23
Resetting the scale	23
Using the external display	24
Set the time.....	24
Alarm function.....	24
Illuminate display.....	25
Recall stored values	26
Change temperature unit	26
Cleaning.....	26
Fault/Remedy.....	27
Technical data.....	27
Declaration of conformity	28
Warranty	29
Disposal.....	29

Scope of delivery

- Bathroom scales
- Extern display
- Battery type AAA, 1.5 V  (4x)

Safety

- This product is not intended to be used by persons (including children aged 8 and above) with reduced or impaired physical, sensory or mental capacities or persons lacking appropriate knowledge or experience unless a person responsible for their safety provides either supervision or appropriate instructions concerning the use of the product before its use.
- Children must not play with the product.
- Cleaning and user maintenance must not be carried out by children unless they are 8 years of age or older and are being supervised.
- Batteries, when swallowed, can be life-threatening. For this reason, the scale and battery should be kept in a place out of reach of children. If a battery has been swallowed, medical assistance must be sought immediately.

DANGER for children

- Keep packaging material away from children. There is a risk of choking if swallowed. Also pay attention that children do not pull the packaging material over their heads.

RISK of fire/burning and/or explosion

- There is a risk of explosion if the batteries are replaced improperly. Replace the batteries only with the same or equivalent battery type. Pay attention to the "Technical data".
- Batteries must not be recharged, reactivated by other means, dismantled, thrown into fire, or short-circuited. Failure to comply with this requirement can lead to explosion and leakage, causing gases to escape!
- Avoid contact with skin, eyes, and mucosa. If contact is made with the battery fluid, rinse the affected areas immediately with plenty of clean water and see a doctor straight away.

WARNING of risk of injury

- If the scale falls down and shatters, there is a risk of injury from glass splinters.
- Do not step onto the outer edge of any one side of the scale. Risk of tipping!
- Do not step onto the scale with wet feet or when its surface is damp. Risk of slipping!

BEWARE of material damage

- Place the scale on a stable, even surface.
- Do not overload the scale. See section "Technical data". Place yourself on the scale exclusively for weighing purposes.
- Protect the scale and the display from strong impacts and vibrations.
- Remove the batteries from the item when they are used up or when you will not be using the item for a longer period. This prevents damage that may occur due to leakage.
- Ensure the polarity (+/-) is correct when inserting and replacing the batteries.
- Do not subject the batteries to any extreme conditions, by – for example – storing them on radiators or under direct sunlight. Increased risk of leakage!
- Clean the battery and device contacts, as needed, before inserting the battery.
- Do not make any modifications to the product. Only allow repairs to be carried out by a specialist workshop. Repairs carried out improperly can considerably endanger the user.
- Do not place the scale in the vicinity of devices that may emit electromagnetic radiation, such as mobile phones, wireless devices, or the like. Doing so may lead to incorrect display indications or malfunctions.
- There are feet on the underside of the scale. Since floor surfaces are made from various materials and are treated with diverse care products, it cannot be completely excluded that some of these substances contain ingredients that may adversely affect and soften the feet of the device. If necessary, place a non-slip pad under the scale.


Intended use

The scope of delivery consists of a personal scale and the external display.

- The personal scale is intended for weighing persons in a private environment. It is not suitable for use in clinics, medical practices, etc. The maximum load capacity of 180 kg must not be exceeded.
- The external display shows the measured weight in conjunction with the scale, whereby the data is transmitted wirelessly.
- Otherwise, the external display serves as a clock with an alarm function and also displays temperature and humidity. The external display is not suitable for commercial use.


Overview

Scale

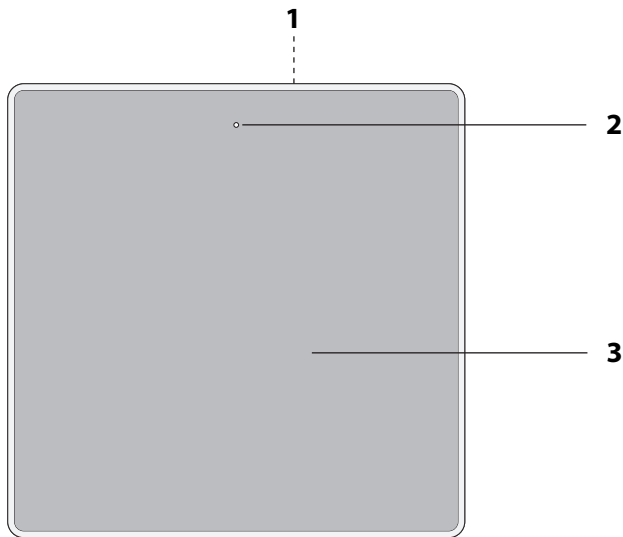
- 1 Battery compartment (bottom)
- 2 Display
- 3 Weighing surface

External Display – Front

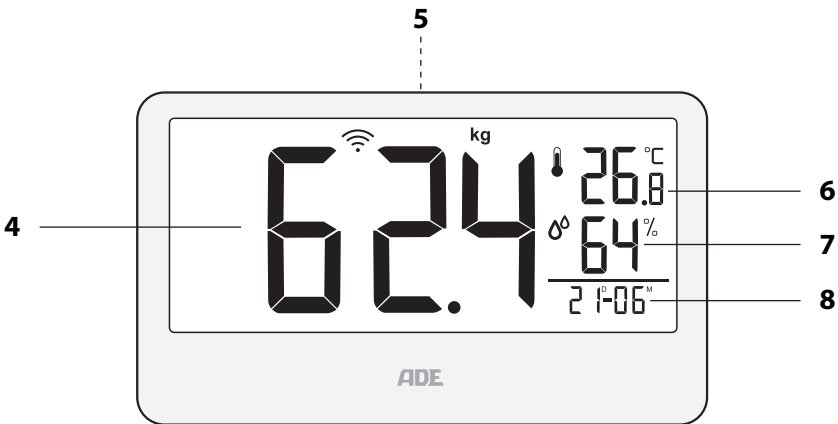External Display – Back

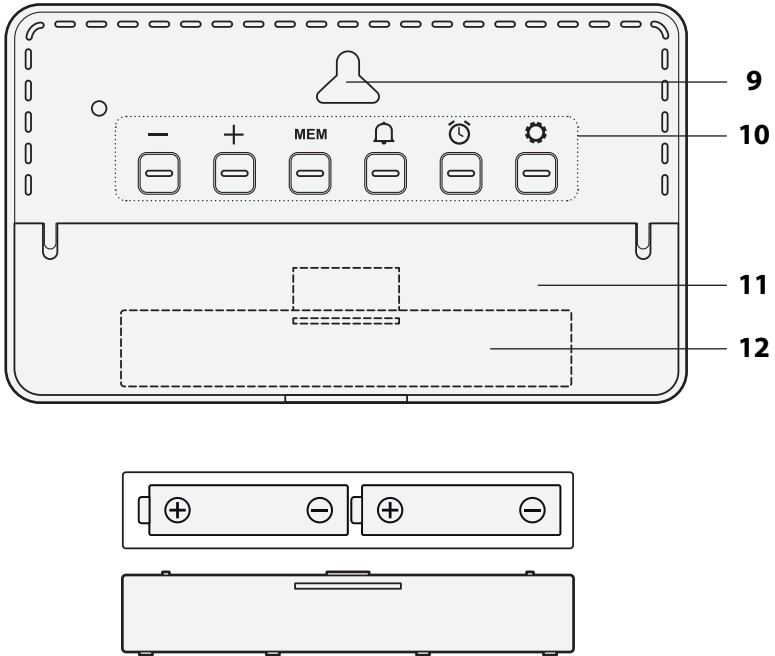- 4 Time or display of measured body weight
- 5 Button ☀ / **ZZ** for snooze function and light
- 6 Temperature in °C or °F
- 7 Relative humidity expressed as a percentage
- 8 Date
- 9 Suspension eyelet
- 10 Setting buttons:
 - + Increase or decrease values, change units of measurement
 - MEM** Retrieve saved values
 - 🔔 Turn the alarm function on and off
 - 🕒 Set alarm times
 - ⚙ Enable settings
- 11 Foldable stand
- 12 Battery compartment

Start-up

Inserting/changing batteries

1. Slide off the battery compartment cover in the direction of the arrow.
2. Insert two batteries as shown on the bottom of the battery compartment. Observe correct polarity (+/-).
3. Slide the battery compartment cover back on.

Insert/replace batteries in the external display

1. Fold up the stand on the back.
2. Remove the battery compartment cover.
3. Insert two batteries as shown on the bottom of the battery compartment. Observe correct polarity (+/-).
4. Replace the battery compartment cover. It must audibly snap into place.

Suggestions for changing the battery:

- A battery change is necessary if  or nothing at all is displayed, see "Malfunction / Remedy".
- Only use the type of batteries specified in "Technical data".
- Always change all batteries at the same time.
- Dispose of the used batteries in an environmentally responsible way (see "Disposal").

Setting up the scale

- Place the scale on an even, firm and dry surface. Do not place the scale on a carpet.

Set up/hang external display

The external display has a hanging eyelet and a fold-out stand on the back for secure placement or wall mounting.

Using the scale

Weighing

WARNING against injuries

- Do not step onto the scale with wet feet or when its surface is damp. Risk of slipping!

1. Hold the external display in your hand or position it so that you can read it well when weighing.
2. Position yourself on the scale. Stand still and do not hold on to anything.
The radio transmission signal lights up and  is displayed as the symbol for data transmission. After that, your measured weight is displayed in the selected unit.
3. To change the unit of measurement, press the button + on the external display (kg > lb > lb st).

A few seconds after you leave the weighing surface, the scale switches off automatically. The measured body weight is displayed for a while on the external display, then the time is shown again.

- Weigh yourself with as little clothing on as possible.
- Try to always weigh yourself at the same time and under consistent conditions in order to obtain comparable values.

Resetting the scale

Reset the scale to its default settings if unusual or obviously incorrect values are shown on the display.

1. Remove the batteries from the battery compartment of the external display and wait a moment.
2. Reinsert the batteries.
Pay attention to the correct polarity of the batteries (+/-).

After switching the scale on, you can weigh yourself as usual.

Using the external display

The external display offers the following functions:

- Display of measured weight, time, temperature, humidity
- Alarm function with snooze
- Display lighting

Set the time

- If more than 20 seconds pass between two button presses, the setting process is automatically canceled and must be repeated.
- If you hold down $-$ or $+$, the displayed values will scroll faster.

1. Hold down the button  until the hour display blinks.
2. Set the current hour using the buttons $-$ / $+$.
3. Press . The minute display blinks.
4. Continue as described above:
 - Press the button $-$ / $+$ to set the desired value.
 - Press  each time to save a setting.
5. Make the following settings:
 - Minute
 - Year
 - Month
 - Day
 - Date format (day/month or month/day)
6. Finally, press  to finish the process.

Alarm function

Set alarm time

1. Press button  once. The alarm time and AL (= ALARM) are displayed.
2. Hold down the button  until the hour digit blinks.
3. Set the desired alarm hour using the buttons $-$ / $+$.
4. Press . The minute display for the alarm time blinks.
5. Set the desired alarm minute using the buttons $-$ / $+$.
6. Press . The set snooze time blinks (default setting: 5 minutes).

7. Set the desired snooze time (time between two alarms) using the buttons $-$ / $+$. If you choose OFF, no automatic alarm repeat occurs. Finally, press  to complete the process.  is displayed, the alarm function is activated.

To later check the set alarm time, briefly press the button  once.

Being woken

At the set time, you will be woken by an alarm signal, which becomes more intense as the wake-up time increases. In addition, the alarm symbol  blinks on the display.

Interrupt alarm for 24 hours

- Press any button (except  / **z^z**) to interrupt the alarm. After 24 hours, the alarm is triggered again.

Auto-stop function

If you do not press a button, the alarm switches off automatically after 2 minutes. After 24 hours, the alarm is triggered again.

Snooze function (SNOOZE)

- Touch  / **z^z** once when the alarm signal sounds.

The alarm signal stops according to the set interval (default setting: 5 minutes).

While the snooze function is active, **z^z** and  blink on the display.

Then the device switches on again and wakes you up again. This process can be repeated several times.

Switch alarm function on/off

From the normal time display:

- Press the button  once.

As soon as the alarm symbol  disappears, the alarm function is deactivated and you will no longer be woken. By pressing  again, the alarm function is reactivated.

Illuminate display

- Briefly tap  / **z^z** to switch on the display light for a few seconds.

Recall stored values

- Repeatedly press **MEM** to see the stored maximum values of the last 24 hours for measured weight, temperature, and humidity (MAX > MIN > currently measured values).
- To delete the stored values, press and hold **MEM** for a few seconds until the display is cleared. Otherwise, the stored data is automatically deleted at midnight every day.

Change temperature unit

From the normal time display:

- Press the button — to switch between the units °C and °F.

Cleaning

NOTE on material damage

- Do not immerse the scale in water or other liquids.
- Never use harsh, scouring or abrasive cleaning agents or the like to clean the scale. Otherwise, the surface could be scratched.
- Wipe the scale as needed with a slightly damp soft cloth.
- Disinfect the weighing surface from time to time using a suitable disinfecting agent (available at specialist retailers and pharmacies).
- Dust the external display if necessary with a soft, dry cloth. The glass can also be cleaned with a commercially available glass cleaner.

Fault/Remedy

Fault	Possible cause
No function	Are the batteries empty or inserted the wrong way around?
The scale appears to be displaying wrong weighing results	- Is the scale located on a firm, even surface? - Is a device that can cause electromagnetic disturbances located in the vicinity of the scale? - Remove the batteries from the scale and wait a moment. Reinsert the battery.
The display shows 'O-LD'.	The maximum capacity of the scale (180 kg) has been exceeded.
The display shows  .	The batteries are exhausted and need to be replaced: - Display  top edge: Replace the batteries in the scale. - Display  next to the humidity display: Replace the batteries in the display.
The temperature reading appears to be too high.	Was the external display exposed to direct sunlight?
The display shows LL.L or HH.H.	The temperature is below or above the measurable range.
The display shows 19% or 96%.	The humidity is below or above the measurable range.

Technical data

Scale

Model:	Digital bathroom scale with external display BE 2500-x "Ronja"
Batteries:	2x AAA / 1,5 V 
Rated current:	25 mA
Measuring range:	max. 180 kg / min. 5 kg
Variance:	1 % +/- 100g

Display scaling:	100 g
Weight:	ca. 1,5 kg
Scale dimensions:	approx. 310 x 310 x 22 mm

External display

Batteries:	2x AAA / 1,5 V 
Transmission frequency:	2,47 GHz
Transmission power:	-39 dBm
Transmission range:	approx. 5 metres
Rated current:	25 mA
Measuring range	
- Temperature:	-9,9 °C to 50 °C (14,2 °F – 122 °F)
- Relative humidity:	20 % – 95 %
External display dimensions:	approx. 143 x 84 x 22 mm
Distributor:	GREND S GmbH
(no service address)	Stahltwiete 23 22761 Hamburg, Germany

Our products are constantly being developed and improved.
Therefore, we reserve the right to make changes to the design and specifications without further notice.

Declaration of conformity

GREND S GmbH hereby declares that the Digital Bathroom Scale BE 2500-x complies with Directives 2014/53/EU and 2011/65/EU.

The declaration will lose its validity if any modification is made to the device without our consent.

Hamburg, july 2025
GREND S GmbH

Warranty

GREND'S GmbH warrants for 2 years from the date of purchase that it will remedy any defects due to faulty material or workmanship free of charge by repair or replacement. In the event of a warranty claim, please return the product with the proof of purchase (specifying the reason for claim) to your dealer.

Disposal

Disposing of the packaging

Dispose of the packaging according to type. Dispose of cardboard with waste paper and of films with recyclable materials.

Disposing of the product

Dispose of the product in accordance with the regulations that apply in your country.

Devices must not be disposed of in normal household waste.

At the end of their life, the products must be disposed of in an appropriate way. In this way, valuable materials contained in the device will be recycled and the burden on the environment avoided. Hand in the old device at a collection point for electronic waste or at a recycling depot.

For more information, contact your local waste disposal firm or your local administrative body.

Recyclable and disposable batteries do not belong in household waste.

All batteries must be handed in at a collection point in your community, suburb, or at the retailer. Such batteries can thus be disposed of in an environmentally friendly way.

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Hierbij verklaart GREND S GmbH dat de digitale weegschaal BE 2500-x conform de richtlijnen 2014/53/EU en 2011/65/EU is.