Bedienungsanleitung
Operating Manual
Manuel d'utilisation
Bedieningshandleiding
Manual de usuario
Manual de instruções

ADE Special

Teleskop- Längenmessstab
Telescopic Height Measuring Rod
Toise de mesure téléscopique
Telescopische lengtemeetstaaf
Estadiómetro telescópico
Estadiómetro telescópico

MZ10023-1
MZ10023-2

MZ10023-261013-Rev008-UM-de-e-fr-nl-es-pt

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Intended use and safety instructions

This symbol means: "Consult operating instructions". To achieve accurate measuring results, read the instructions for use carefully before use and follow all instructions provided therein. Keep the instructions for use for future reference and in reach for the user.

Your ADE measuring rod is a high-quality product. To achieve the best possible results, please carefully follow the instructions in this user manual.

The measuring rod is designed to measure the height of adults and children. It can be used either by mounting it on the wall or in combination with an ADE column scale.


Safety instructions

- Follow all instructions and notes in this operating manual.
- Before assembling and using the measuring rod, check it for damage. If damaged, do not use the device due to the risk of injury.
- When the measuring rod is not being used for measuring, always swivel the measuring tongue down completely to prevent any risk of injury.
- Only attach the measuring rod to ADE column scales or to a wall using a measuring stick and spirit level.
- In the installation area, the wall must be level and plumb.
- In the installation area, the floor on which the person to be measured will stand must be level and horizontal.
- Skirting boards or plinths must not protrude further into the room than the installed height measuring rod.
- If the ceiling height is too low, the measuring range may be restricted.

Contents of the package / Fitting materials

Carton content

- Telescopic measuring rod
- Mounting material for the two possible installation types (see image)
- User manual

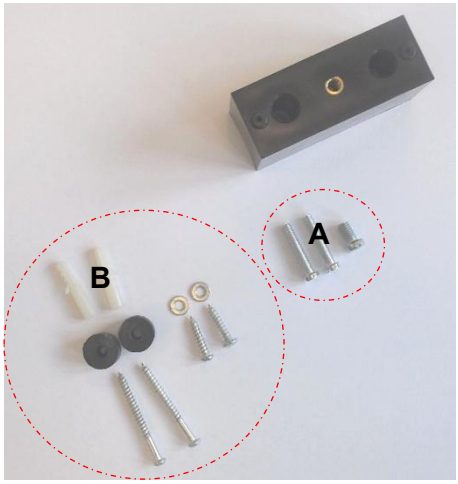Spacer block for mounting to the front or rear side of the column of a column type scale.

A: Fixation screws for fixing the measuring rod to the scale column

B: Screw set for wall mounting incl. wall dowels, screws and washers

A: Column scale assembly	B: Wall assembly
• 1 spacer block for mounting to the scale column • 1 M6 screw • 2 M5 screws	• 2 spacer washers (black plastic) • 2 screws, short (for mounting to wood or metal walls) • 2 dowels S6 • 2 long screws (for mounting to concrete walls) • 2 washers


Mounting

Mounting the Measuring Rod to the Scale Column

Your ADE column scale is equipped with two threaded holes on the back of the column. These allow mounting of the spacer block for the measuring rod.

Position the spacer block over the threaded holes on the column and fasten it into the fixing bores **1** using the two M5 screws. Tighten both screws firmly.

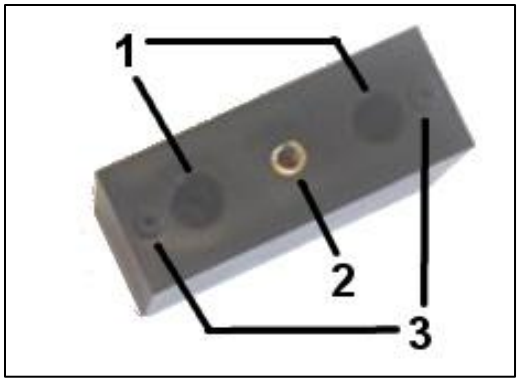1: fixing bores

2: thread for fixing screw

3: locking pins


The spacer block additionally includes a threaded hole **2** for the measuring rod's mounting screw, as well as two positioning pins **3** spaced symmetrically. At the lower end of the measuring rod, on the printed side (measuring scale), there are corresponding holes for alignment.


Position the measuring rod onto the spacer block so that the locking pins align with the corresponding holes. Secure the measuring rod using the M6 screw in the lower mounting hole. Tighten the screw firmly.

Mounting the Measuring Rod to the Wall

Tools needed: folding rule, pencil, level, drill, 6 mm drill bit, screwdriver

1. Measure 56 cm from the floor (or from the surface where the person will stand later) on the wall, accurate to the millimeter, and mark the spot. Drill the first hole there and insert the dowel.
2. Extend the rod all the way. Insert the first screw with a washer from the front through the topmost of the four holes at the lower end (directly below the start of the scale at 60 cm). Place the black spacer washer on from the back and lightly screw the screw into the wall plug (the rod must still be able to rotate).
3. Hold a level against the side of the rod and align the rod vertically. Mark the center of the hole on the wall through the hole (at the 122.5 cm mark on the scale).
4. Swivel the rod to the side or remove it briefly. Drill the second hole and insert the dowel.
5. Swivel the rod back into place. Install the second screw (with the washer in front and the spacer disc in back).
6. Tighten both screws securely.

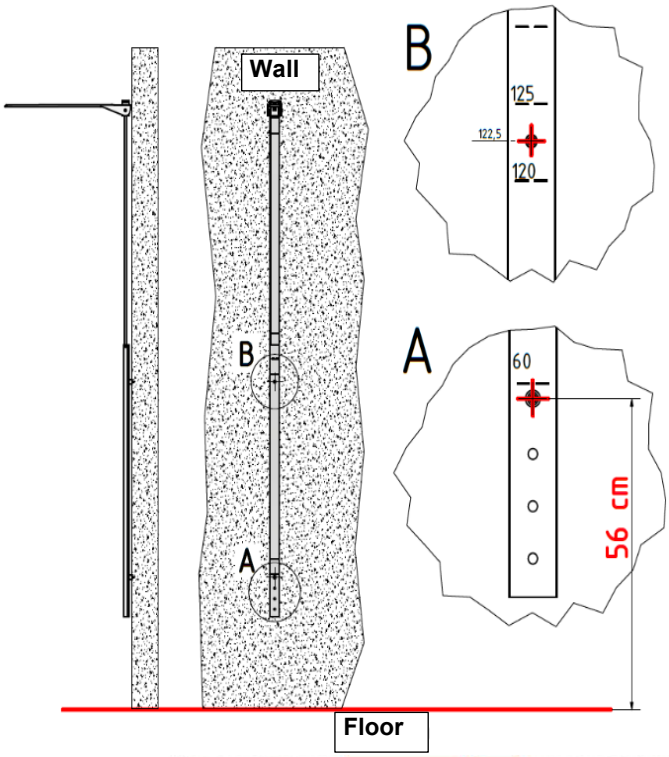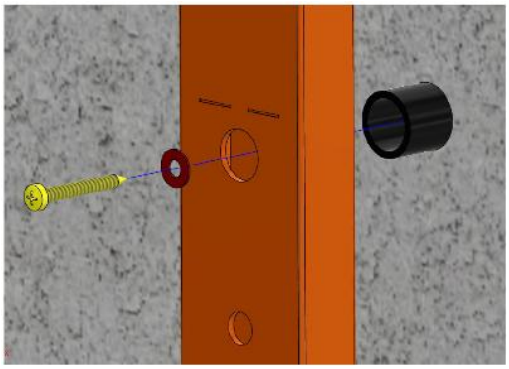Operating


1. Move the measuring tongue to the horizontal position. It will audibly and noticeably lock into place.

2. Raise the measuring tongue above the person to be measured.
3. Have the person stand directly under the measuring tongue.
4. Lower the measuring tongue until it rests flat on the person's head.
 (illustration on principle; ill. showing different ADE product)

5. Read the measurement at the red line located between the lower and upper part of the measuring rod. For measuring persons under 134 cm press the black button in the center of the measuring tongue and slide the tongue further down onto the lower section of the measuring rod. Read the measurement at the red line located at the lower edge of the measuring tongue.

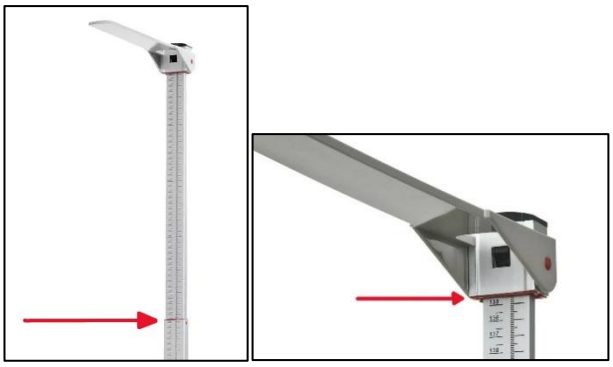6. When not in use, swivel the measuring tongue downward to prevent injury. A slight resistance is noticeable while swiveling the tongue.

Cleaning and Maintenance

For cleaning the measuring rod use a damp, soft and lint-free cloth. If necessary, use mild detergents only (e.g. household cleaner).

Do not use any aggressive, abrasive, acid or chlorine containing detergents or disinfectants.

If a surface disinfection becomes necessary, use an alcohol-based disinfectant and follow all dosing and handling instructions given by the respective manufacturer of the agent.

Warranty – Liability

ADE will repair or replace this product within a period of 2 years from date of delivery applying to defects occurred due to poor material or workmanship (presentation of invoice is necessary). Warranty does not include natural wear, damages caused by incorrect or negligent handling, overuse, unsuitable or improper storage or inappropriate transport methods.

Customers abroad are kindly requested to contact their local dealer for any warranty claim, customer service and spare parts.

Disposal of the equipment:

If the measuring rod has reached its end of service life, dispose of it properly and professionally in accordance with the country-specific laws and local regulations.

Technical Data

Measuring Range:	600 – 2100 mm
Graduation:	1 mm
Dimensions (WxHxD):	55 x 940 x 34 mm
Weight:	approx. 700 g

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ADE Germany GmbH